


PINNA[®]

Your Swiss Option

VOLUTUS

INVADUR[®]

VERTEBRAL
BODY
REPLACEMENT

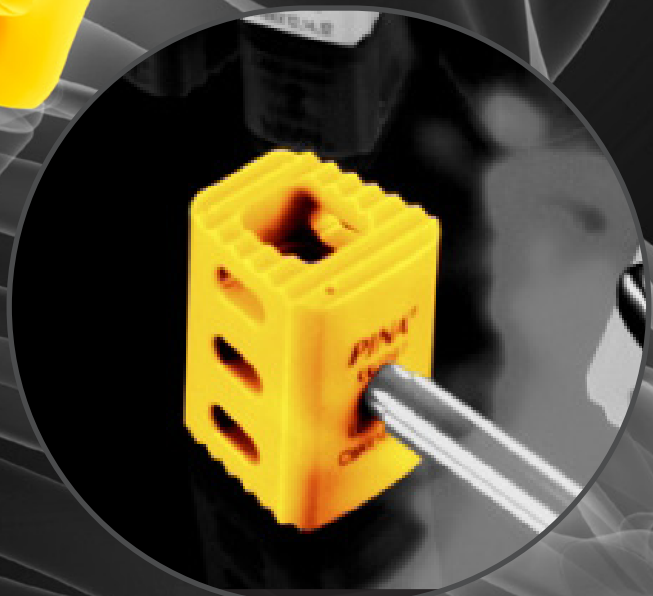



VOLUTUS INVADUR[®] FEATURES

The PINA[®] VOLUTUS Invadur[®] is indicated for the replacement of cervical vertebral bodies in skeletally mature patients following uncuser-preserving corpectomy at a level from C3 to C7. The implants can be implanted through a lateral anterior approach to the cervical spine. Their shape is adapted to the morphology of the vertebral bodies to be treated and the prescribed surgical techniques. They have a cavity to accommodate autogenous bone material. The vertebral body replacement is available in various designs and sizes. The VOLUTUS Invadur[®] is made of polymers from the PAEK group. This specially made Invadur[®] material is X-ray transparent, bio-, MRI-, and CT-compatible.

FEATURES

- Implants can be combined using the adapter inserted into the central canal
- The lordosis angle is 4,1° or 3,6° - a plane-parallel version is also available
- Eight different heights (12 mm - 26 mm) with an increment of 2 mm
- In case the implants are combined, the lordosis angle is 8,2° or 7,2°
- The easy-to-separate trial set makes it easier to determine the size
- Three-dot x-ray marker for visualization during imaging




VOLUTUS INVADUR[®] IMPLANTS

UNSTERILE

ARTICLE NUMBER	ARTICLE DESCRIPTION	SIZE (WIDTH X LENGTH)	HEIGHT	ANGLE
CWKK 12.14.XX	Volutus Invadur	12 x 14	XX	4,1°
CWKK 12.14.14 P	Volutus Invadur planparallel	12 x 14	14	0°
AK 12.14.	Stacking adapter	12 x 14	-	-
CWKK 14.16.XX	Volutus Invadur	14 x 16	XX	3,6°
CWKK 14.16.14 P	Volutus Invadur planparallel	14 x 16	14	0°
AK 14.16.	Stacking adapter	14 x 16	-	-

STERILE

ARTICLE NUMBER	ARTICLE DESCRIPTION	SIZE (WIDTH X LENGTH)	HEIGHT	ANGLE
CWKPOS 12.14.XX	Volutus Invadur	12 x 14	XX	4,1°
CWKPOS 12.14.14 P	Volutus Invadur planparallel	12 x 14	14	0°
AKPOS 12.14.	Stacking adapter	12 x 14	-	-
CWKPOS 14.16.XX	Volutus Invadur	14 x 16	XX	3,6°
CWKPOS 14.16.14 P	Volutus Invadur planparallel	14 x 16	14	0°
AKPOS 14.16.	Stacking adapter	14 x 16	-	-

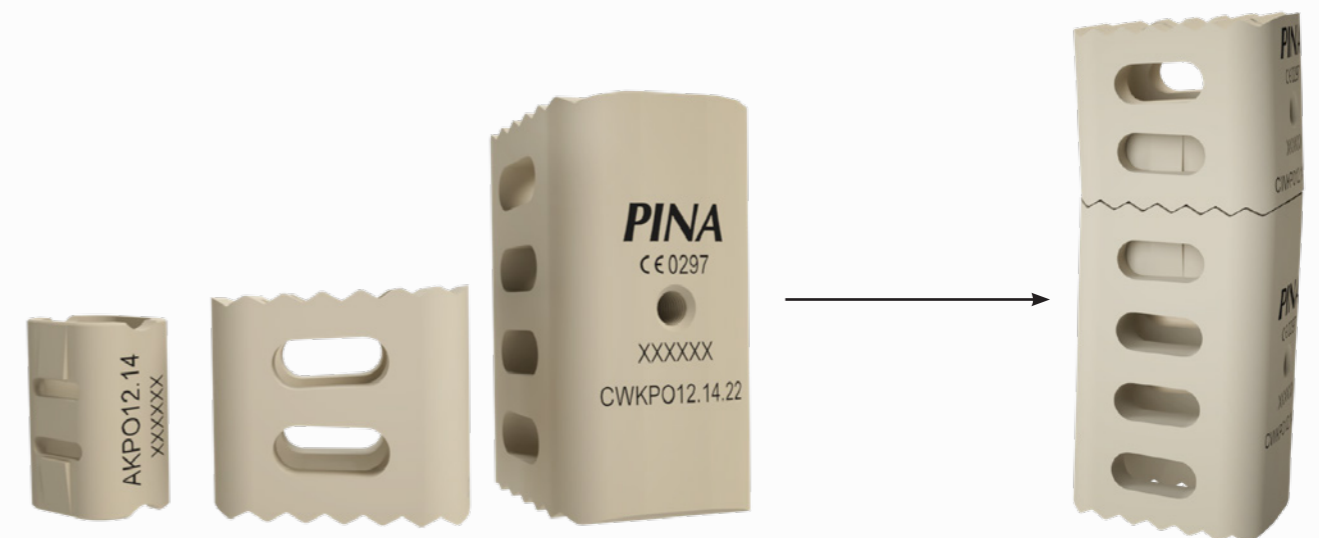



VOLUTUS INVADUR®

MORE INFORMATION

INSTRUMENTS

- Clearly arranged instruments
- Blue colored trial implants



Further information at
info@bricon.com
and
+49 (0) 7461 9336-0

ADDITIONAL SYSTEMS

Other systems for the cervical and lumbar spine can be found under



bricon.com/en/produkte-en/wirbelsaule-en/

PINA®

Your Swiss Option

PINA Medizintechnik Vertriebs AG
Neuwiesenstrasse 15
8400 Winterthur / Switzerland

Tel.: +41 (44) 586 52 - 62

info@pina-med.ch
www.pina-med.ch

BRICON

Distributed by:

Bricon GmbH
Eisenbahnstrasse 100
78573 Wurmlingen / Germany

Tel.: +49 (0) 7461 9336 - 0

info@bricon.com
www.bricon.com